



Fracture Clinic Product Catalogue

KSA

Refining our product range for over thirty years
we proudly introduce our catalogue specifically
for clinicians in trauma or fracture clinics.

Please take a look at the products we offer,
feel free to request any samples for a
closer look at items of interest.

Contents



Product

OptimaCast	2
OptimaCast Soft	3
Optima Splint	4
Optima Pad & Stockinette	5
Scissors	6
Cast Removal	7
Hook & Loop Fasteners	8-9
Padding	10
Immobilisers	11
Footwear	12-13
Wrist Braces	14-15
Fracture Braces	16

OptimaCast

- OptimaCast and OptimaCast Poly is used in the construction of most common orthopaedic casts.
- Consisting of knitted fibreglass or polyester, OptimaCast provides a rigid cast which will immobilise any required area.
- Both impregnated with a water-activated polyurethane resin.
- OptimaCast Poly is lightweight, strong and radioluscent.
- Compared with its fibreglass counterpart OptimaCast Polyester is more lightweight and can be completely burnt and therefore more environmentally friendly.



1. Wear protective gloves.
2. Apply stockinette over the desired area to be cast.



3. Immerse OptimaCast in water for approx 5 seconds.
4. Squeeze the tape tightly to remove excess water.



5. Wrap the desired area with the tape overlapping itself by one third to one half each time. Do not apply too tightly.



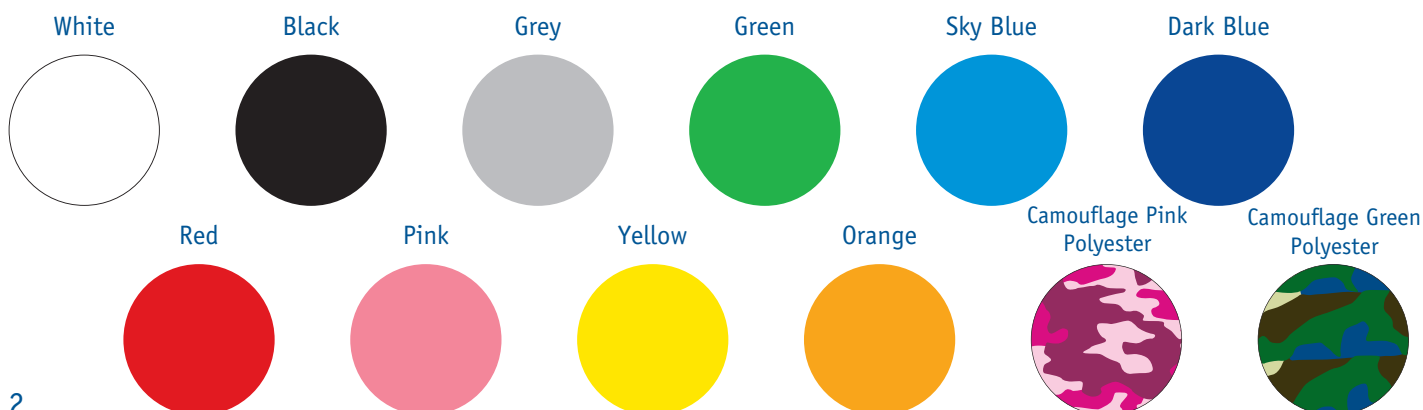
Moulding and forming should be done at this time.

Setting time is approximately 3-5 minutes.

Functional strength will be achieved in 20 minutes.

Code	Description
CP001	OptimaCast Polyester 1" (2.5cm) x 1.8m
CP002	OptimaCast Polyester 2" (5cm) x 3.6m
CP003	OptimaCast Polyester 3" (7.5cm) x 3.6m
CP004	OptimaCast Polyester 4" (10cm) x 3.6m
CP005	OptimaCast Polyester 5" (12.5cm) x 3.6m

Code	Description
CF001	OptimaCast Fibreglass 1" (2.5cm) x 1.8m
CF002	OptimaCast Fibreglass 2" (5cm) x 3.6m
CF003	OptimaCast Fibreglass 3" (7.5cm) x 3.6m
CF004	OptimaCast Fibreglass 4" (10cm) x 3.6m
CF005	OptimaCast Fibreglass 5" (12.5cm) x 3.6m



OptimaCast Soft

- A knitted polyester based fabric impregnated with a polyurethane resin.
- On exposure to moisture a chemical reaction will make the tape semi-rigid.
- The splint/cast will then retain it's shape yet allow very small movement aiding the rehabilitation process.

FLEXIBLE

- OptimaCast Soft is intended for use when support is required. It is not for immobilisation alone.
- Ideal for paediatric corrective serial casting, secondary casting for various orthopaedic problems.

EASILY REMOVED

- OptimaCast Soft can be unwrapped and removed with bandage scissors.



1. Wear protective gloves.
2. Apply stockinette over the desired area to be cast.



3. Immerse OptimaCast Soft in water for approx 5 seconds.
4. Squeeze the tape lightly to remove excess water.




5. Wrap the desired area with the tape overlapping itself by one third to one half each time. Do not apply too tightly (OptimaCast Soft requires 3-4 layers to achieve the desired level of immobilisation).



Moulding and forming should be done at this time.
Setting time is approximately 3-5 minutes.
Functional strength will be achieved in 20 minutes.

Code	Description
CS001	OptimaCast Soft 1" (2.5cm) x 1.8m
CS002	OptimaCast Soft 2" (5cm) x 3.6m
CS003	OptimaCast Soft 3" (7.5cm) x 3.6m
CS004	OptimaCast Soft 4" (10cm) x 3.6m
CS005	OptimaCast Soft 5" (12.5cm) x 3.6m

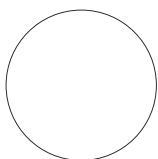White




Black




Grey



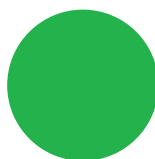
Green



Sky Blue

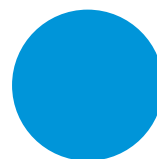


Dark Blue

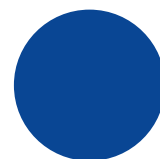


Red



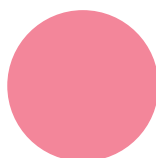Pink



Yellow

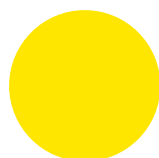

Orange



Optima Splint

- Fibreglass material encased in soft anti-bacterial cloth. Used for protecting broken skin and immobilising wounded region.
- Available in pre-cut and roll form.



1. Select the suitable width and length according to the region to be treated.
2. Soak the splint for 5-6 seconds in water and remove.
3. Roll the splint and squeeze 2-3 times to remove moisture.
4. The splint can then be rolled in a towel to absorb excess water.
5. Mould this to the desired area and wrap in an elasticated bandage.

Optimasplint will activate on contact with moisture in the air. Ensure the splint rolls are properly closed after use.

OPTIMA SPLINT PRECUT TYPE

Code	Description
SF210	Optima Splint 2" (5cm) x 10"
SF312	Optima Splint 3" (7.5cm) x 12"
SF335	Optima Splint 3" (7.5cm) x 35"
SF415	Optima Splint 4" (10cm) x 15"
SF430	Optima Splint 4" (10cm) x 30"
SF530	Optima Splint 5" (12.5cm) x 30"
SF545	Optima Splint 5" (12.5cm) x 45"
SF630	Optima Splint 6" (15cm) x 30"
SF645	Optima Splint 6" (15cm) x 45"

OPTIMA SPLINT ROLL TYPE

Code	Description
SR001	Optima Splint Roll 1" (2.5cm)
SR002	Optima Splint Roll 2" (5cm)
SR003	Optima Splint Roll 3" (7.5cm)
SR004	Optima Splint Roll 4" (10cm)
SR005	Optima Splint Roll 5" (12.5cm)
SR006	Optima Splint Roll 6" (15cm)

Optima Pad & Stockinette

- Undercast padding for cushioning and protection under synthetic casts and plaster of paris bandage.
- Made from 100% finest quality polyester.
- Keeps the skin dry and comfortable for the patient during long term wear.
- Smooth surface and good air permeability.



OPTIMA PAD

Code	Description
ROCPad2	2" (5cm) x 360cm Pad
ROCPad3	3" (7.5cm) x 360cm Pad
ROCPad4	4" (10cm) x 360cm Pad
ROCPad6	6" (15cm) x 360cm Pad

- Made from skin friendly, finest quality material.
- Use under synthetic cast or plaster of paris bandage to protect against irritation during the long casting process.
- Packed and sold in easy to use dispenser cartons.

STOCKINETTE

Code	Description
OPNet002	2" (5cm) x 23m Stockinette
OPNet003	3" (7.5cm) x 23m Stockinette
OPNet004	4" (10cm) x 23m Stockinette
OPNet006	6" (15cm) x 23m Stockinette



Scissors

BANDAGE SHEARS

With slanted, stainless, serrated blades with blunt tips and needle destroyer. Plastic handles. Available with nonstick coating. Autoclavable.

Code	
RA73386	Bandage Shears



LISTER BANDAGE SCISSORS

Made from steel available in 5^{1/2}" and 8". Large finger ring.


Code	
RA98SCS1-26	Scissors, Large Finger Ring, 5 ^{1/2} "
RA98SCS1-19	Scissors, Large Finger Ring, 8"



WIRE CUTTING SCISSORS

Nibbles Angled with One Serrated Blade, 4^{3/4}".

Code	
RA98SCS71-1	Scissors, wire cutting, angled, economy

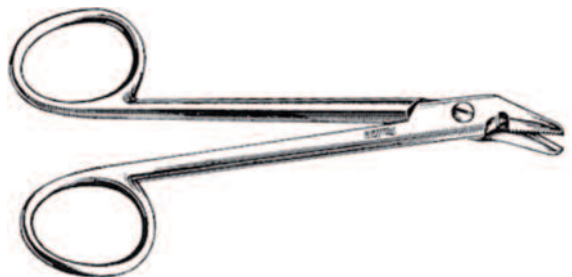

SCISSOR HOLSTER

A solution to having too many scissors and not enough pockets! The scissor holster has two wide slide in pockets and three slide in pockets. All beneath a fold over cover. Attachable to belts up to 2^{1/2}" wide.

Code	
RA932401	Scissor Holster



Cast Removal

BEESON CAST SPREADER

1½" Wide blades, 12".

Code	
RA98CA5S5-1	Economy Beeson Cast Spreader, 12"
RA98CAS5-2	Miltex 27-3090 Beeson Cast Spreader, 12"



STRYKER TYPE CAST SPREADER

3 prongs, spring action, Miltex, 9".

Code	
RA98CAS8-2	Miltex 27-3100 Stryker Type Cast Spreader, 9"



WOLF-BOEHLER CAST BREAKER

Deep serrations, 15mm wide blades.

Code	
RA98CAS9-1	Wolf-Boehler Cast Breaker, Economy, 10"
RA98CAS9-2	Wolf-Boehler Cast Breaker, Economy, 7"
RA98CAS9-3	Wolf-Boehler Cast Breaker, Miltex, 10"





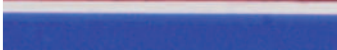

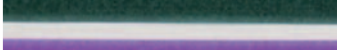


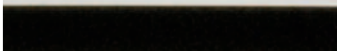
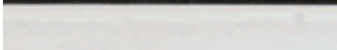












HENNING CAST SPREADERS

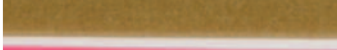2" wide blades, 11".

Code	
RA98CAS6-1	Economy Henning Cast Spreader, 11"
RA98CAS6-3	Miltex Henning Cast Spreader, 11"



Hook & Loop Fasteners

	DARK GREY - 33
	BURGUNDY - 168
	ROYAL BLUE - 55
	KELLY GREEN - 65
	LILAC - 298
	JADE - 659
	BLACK - 90
	WHITE - 10
	RED - 52
	BEIGE - 20
	FUCHSIA - 171
	NAVY BLUE - 18
	YELLOW - 75
	LIGHT GREY - 32
	BRIT MILITARY GREEN - 07
	FLO. ORANGE - 390
	LIGHT BLUE - 27
	FLO. YELLOW - 123
	BROWN - 25
	FLO. PINK - 155
	FLO. GREEN 231

STANDARD HOOK OR LOOP Supplied in 25 metre rolls


When ordering, please specify the colour, width and whether hook or loop is required.

Width	
10mm (available in white, beige and black only)	£5.40
16mm	£5.70
20mm	£6.30
25mm	£6.90
30mm	£8.10
38mm	£11.40
50mm	£12.00
100mm (available in white, beige and black only)	£30.00



VELSTRAPS® Supplied in Packs of 100

Strap (Length x Width)	Code	
	(Hook + Loop = Length x Width)	
(16") 41cm x 20mm	100H/325L/20VND	£48.90
(20") 51cm x 20mm	100H/435L/20VND	£49.50
(24") 61cm x 20mm	100H/535L/20VND	£52.50

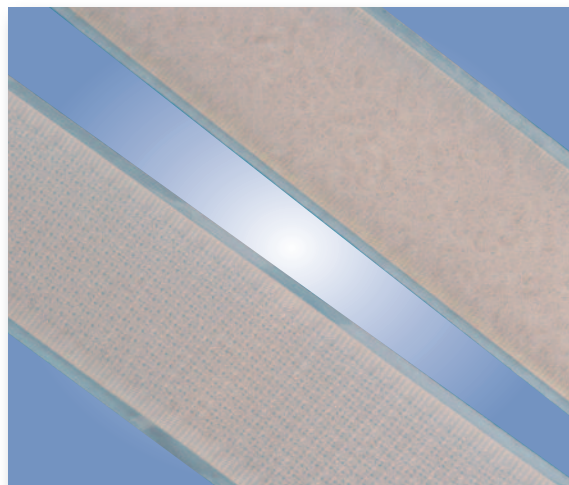

SELF ADHESIVE HOOK OR LOOP Supplied in 25 metre rolls

When ordering please clarify whether hook or loop is required.


Width	Colour	
16mm	White only	£11.10
20mm	White only	£13.60
25mm	White or Black	£13.80
30mm	White or Black	£16.80
50mm	White or Black	£26.00

Hook & Loop Fasteners

DOUBLE SIDED HOOKS

1" x 1½" (2.5cm x 3.8cm)
Double sided Hooks 40 pack

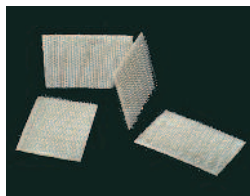RA4240 1" x 1½" (25mm x 38mm)



VELCRO® BRAND SELF ADHESIVE STRAPS

For extra-strong adhesion, warm the hook backing with a heat gun.
Trim loop to length. (white or beige)

RA7120 White 1" x 14" (25mm x 35.5cm)
10 pack

RA7095 Beige 1" x 14" (25mm x 35.5cm)
10 pack



BETA PILE II™ SELF-ADHESIVE STRAPS

Combine stretch-resistant, soft padded Beta Pile with 3" (75mm) of self-adhesive hook. Trim to length or width.
Won't stretch. Beige. 10 pack

RA4464 Beige 1"x16" (25mm x 40.4cm)
RA4465 Beige 2"x16" (50mm x 40.5cm)



SUPER STRAPS®

Super Straps® are tough hook touch and close bonded to soft foam. 1" (25mm) self-adhesive loop on each strap end. 10 pack

RA7652 Super Straps 1"x14" (25mm x 35.5cm)
RA7653 Super Straps 2"x14" (50mm x 35.5cm)



BETA PILE II™ D-RING STRAPS

Soft Beta Pile II combined with 2" (50mm) of hook and loop closure. Nylon D-ring.
Comfortable strap for short-term wear. Beige. 10 pack

RA7100 Beta Pile II D-Ring Straps 2"x14" (50mm x 35.5cm)



SOFT-TOUCH™ D-RING

Soft-Touch™ D-Ring padded. Soft foam between a loop and soft tricot.
No exposed hook to catch.
Cannot be trimmed. Beige. 10 pack



RA7110 Soft-Touch D-Rings, Beige, 1"x10" (25mm x 30.5cm)
RA7111 Soft-Touch D-Rings, Beige, 1"x12" (25mm x 33cm)
RA7112 Soft-Touch D-Rings, Beige, 1"x14" (25mm x 35.5cm)
RA7113 Soft-Touch D-Rings, Beige, 2"x14" (50mm x 35.5cm)
RA7114 Soft-Touch D-Rings, Beige, 2"x18" (50mm x 45.5cm)

D-RING STRAPS

When attached, there will be no exposed hooks to catch on clothing.
Cannot be trimmed. White. 10 pack



RA7055 D-Ring Straps, 1"x10" (25mm x 25cm)
RA7030 D-Ring Straps, 1"x14" (25mm x 35.5cm)
RA7031 D-Ring Straps, 1"x18" (25mm x 45.5cm)
RA7032 D-Ring Straps, 2"x14" (50mm x 35.5cm)

BIG D-RINGS

Large opening. Easy thumb access.
Zytel® nylon. 100 pack

RA4945 Big D-Rings 1"x1½" (25mm x 32mm)



ZYTEL® D-RINGS

Used for fabricating custom D-ring straps.
Nylon. 100 pack

RA4241 Zytel D-Rings, 1" (25mm)
RA4242 Zytel D-Rings, 2" (50mm)



D-RINGS WITH CHAFES

All nylon, easily rivet to splints. 25 pack

RA70833 D-Rings, 1" (25mm)
RA70834 D-Rings, 1½" (32mm)
RA70835 D-Rings, 2" (50mm)



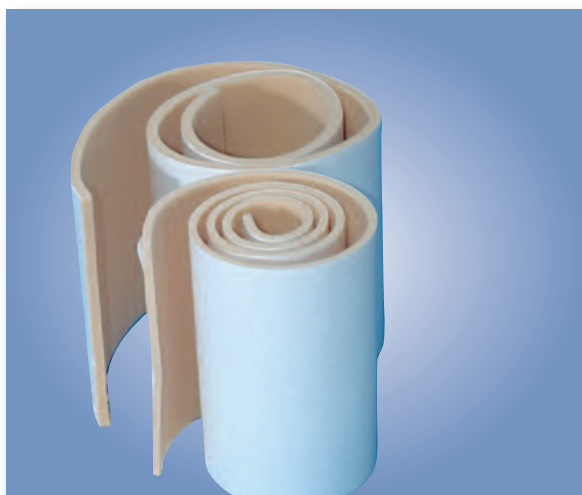
STRAP PADS

Provide a large, cushioned, pressure bearing area. Slots fit straps to 1½" (32mm)W.
Made of Aliplast 4E. 20 apck



RA6422 Strap Pads

Padding



PLASTAZOTE® SELF-ADHESIVE PADDING


Plastazote® self-adhesive padding is made from from Plastazote #1. Can be applied to your splinting material before forming. Simply adhere to the sheet or splint blank before heating. Closed-cell foam. Pink. 6" x 36".

Code	
RA4430	Plastazote Self-Stick, 1/16"x6"x36" (2mmx15cmx91cm)
RA4431	Plastazote Self-Stick, 1/8"x6"x36" (5mmx15cmx91cm)
RA4432	Plastazote Self-Stick, 3/16"x6"x36" (8mmx15cmx91cm)
RA4433	Plastazote Self-Stick, 1/4"x6"x36" (6mmx15cmx91cm)

FELT: PLAIN AND ADHESIVE

Orthopedic-grade felt for lining and padding. Cotton-wool blend for resilience and next-to-skin comfort. Sanitary. Plain felt sold in 10" x 36" rolls. Adhesive backed felt sold in 6" x 2½ yard rolls. White.

Code	
RA4136	Felt, plain, 1/8", 21"x36" roll
RA4137	Felt, plain, 1/4", 21"x36" roll
RA4138	Felt, plain, 1/2", 21"x36" roll
RA4139	Felt, adhesive-backed, 1/8", 6"x2.5yd roll
RA4140	Felt, adhesive-backed, 1/4", 6"x2.5yd roll

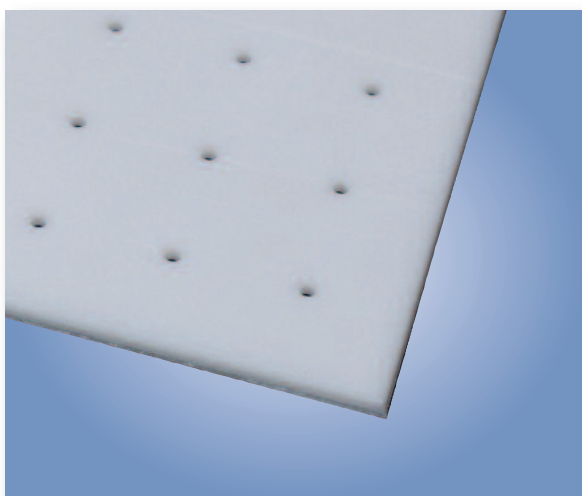

MALLUSPLINT

Perfectly suited to fabricated splints. Malusplint is an easy to use, cost-effective thermoplastic. Simply heat in a water bath until it becomes transparent, then mould to the desired area. Available in a range of colours and % of perforations. Please contact our office for further details.



KSA SHOULDER IMMOBILISER

The shoulder immobiliser is a simple way to relax the arm and shoulder after trauma. Hook and loop fastener with plastic D rings for tightening. Can be used on either left or right arm.

Code	Dimensions
ASHOULD	One Size Fits All



KSA RADCLIFFE KNEE IMMOBILISER

A very popular inexpensive soft brush nylon, adjustable knee immobiliser with 3 panels and 4 lightweight steels to keep the knee in a strong supported position.

Code	Size	
LRAD16	Paediatric	16" (40.5cm)
LRAD18	Youth	18" (45.5cm)
LRAD20	Adult Standard	20" (50.5cm)
LRAD22	Adult Large	22" (55.5cm)
LRAD24	Adult Extra Large	24" (60.5cm)

Footwear



KSA POST-OP SHOE

Made from a soft brush nylon upper and has 8mm micro sole unit. The shoe is supplied with two pieces of hook fastening. Available in three sizes.

Code	Size	Sole Length
FPOSTS	Small	23cm
FPOSTM	Medium	27cm
FPOSTL	Large	30cm



CLASSIC POST-OP SHOE

Durable upper with soft inner lining allows foot to breathe while moulded heel collar contours around heel for better fit. D-ring closure ensures even compression and adjusts to accommodate bulky dressings and edema, while open-toe design helps prevent undue friction on toes. Fits left or right.



Code	Size
RA65466MS	Mens, Small 5-7
RA65466MM	Mens, Medium 7½-9
RA65466ML	Mens, Large 9½-11
RA65466MXL	Mens, X-Large 11½-13
RA65466M2XL	Mens, 2X-Large 13½-15½
RA65467WS	Womens, Small 3-5
RA65467WM	Womens, Medium 5½-7
RA65467WL	Womens, Large 7½-10
RA65468P	Pediatric 1-10

CANVAS ROCKER CAST SHOE

Popular style canvas duck upper. Rocker sole helps promote a natural gait while the open toe and heel allows for more versatility in sizing. Elastic strap around the ankle reduces slippage. Black.

Code	Size	Sole Length
RA65480S	Black, Small	9½" (18cm)
RA65480S/M	Black, Small/Medium	10½" (26.5cm)
RA65480M	Black, Medium	11" (27.5cm)
RA65480M/L	Black, Medium/Large	12½" (31cm)
RA65480L	Black, Large	12½" (32cm)
RA65480XL	Black, Extra Large	13½" (34cm)
RA64440	Teddy Bears	5½" (13cm)
RA64441	Teddy Bears	7½" (19cm)



DARCO® SHOE

It's soft, padded mesh upper with open toe is lightweight and breathable. Adjustable closures accommodate bulky dressings and provide compression. Fits left or right foot.



Code	Size
RA6931	Small Men: 5-7
RA6932	Medium Men: 7½-9
RA6933	Large Men: 9½-11
RA6934	X Large Men: 11+
RA6951	Small Women: 4
RA6952	Medium Women: 4½-6
RA6953	Large Women: 6½-8

Footwear



DARCO® SLIMLINE CAST BOOT (Adult and Child)

Has a stiff sole. May be used as a post-op. Long wearing polyester and EVA upper. Rocker design sole with a ripple finish. Velcro® closures. Blue. Sold individually.

Adult Code	Size	Men's	Women's
RA64716S	Small	up to 4	up to 6
RA64716M	Medium	5-7	6½-8
RA64716L	Large	7½-10	8½+

Paediatric Code	Size	Sole Length
M40DRCBPS	Small	14cm
M40DRCBPM	Medium	17cm
M40DRCBPL	Large	20cm

DARCO® ORTHOWEDGE™

Protects the forefoot by removing pressure from the metatarsal head. Fits left or right foot.

Code	Size	Men's	Women's
RA6090	X-Small	-	up to 4
RA6091	Small	up to 5	4-5½
RA6092	Medium	5-7½	6-8
RA6093	Large	8-10	8½+
RA6094	X Large	11+	-



DARCO® WALKER FX PRO™

Darco® Walker FX Pro™ Ultra low profile design eliminates height difference between walker and normal shoe without altering gait. Padded liner is vented promoting air circulation, keeping the lower leg cool and comfortable. Anatomically contoured "spatula" stirrups comfortably immobilize the ankle while the pressure-moldable, conforming insole absorbs shock. Available in standard and low heights.

Code	Product	Mens	Womens
RA64580	Darco FX Pro Low XS	up to 4	up to 4
RA64581	Darco FX Pro Low Sm	4½-6½	4½-7
RA64582	Darco FX Pro Low Med	7-9½	7½-9½
RA64583	Darco FX Pro Low Lg	10-12	10+
RA64584	Darco FX Pro Low XL	13+	-
RA64585	Darco FX Pro Standard XS	up to 4	up to 4
RA64586	Darco FX Pro Standard Sm	4½-6½	4½-7
RA64587	Darco FX Pro Standard Med	7-9½	7½-9½
RA64588	Darco FX Pro Standard Lg	10-12	10+
RA64589	Darco FX Pro Standard XL	13+	-

Qualcare™ ANKLE BRACE WITH GEL/AIR LINER

Anatomically shaped shell attaches and removes in seconds with included hook and loop straps. Fits nicely in sport or dress shoes. Gel/air combination liner sealed in vinyl.

Code	Size
RA62849	One Size Fits All





KSA NEOPRENE WRIST BRACE

Made from neoprene. has removable aluminium stay held in place by a matching neoprene cover.
Hook and loop fastening. Available in blue/navy or beige.

Size	Beige left	Beige right	Blue left	Blue right	Wrist measurement
Extra small	WRNEOBEXSL	WRNEOBEXSR	WRNEOBLXSL	WRNEOBLXSR	14cm
Small	WRNEOBESL	WRNEOBESR	WRNEOBLSL	WRNEOBLSR	16cm
Medium	WRNEOBEML	WRNEOBEMR	WRNEOBLML	WRNEOBLMR	18cm
Large	WRNEOBELL	WRNEOBELR	WRNEOBLLL	WRNEOBLLR	19cm
X Large	WRNEOBEXLL	WRNEOBEXLR	WRNEOBLXLL	WRNEOBLXLR	20cm

Note: Also available in sizes 2XL and 3XL



KSA NEOPRENE WRIST BRACE WITH THUMB

This wrist and thumb brace is designed to act as a Thumb Spica and a wrist brace all in one. Providing the maximum amount of support possible to the hand, wrist and thumb at the same time. The design of the wrist brace allows the thumb to be held in the pen grip position. Extensively used in A & E Departments for suspected Scaphoid injuries.

Size	Beige left	Beige right	Blue left	Blue right	Wrist measurement
Extra small	WRNTBEXSL	WRNTBEXSR	WRNTBLXSL	WRNTBLXSR	14cm
Small	WRNTBESL	WRNTBESR	WRNTBLSL	WRNTBLSR	16cm
Medium	WRNTBEML	WRNTBEMR	WRNTBLML	WRNTBLMR	18cm
Large	WRNTBELL	WRNTBELR	WRNTBLLL	WRNTBLLR	19cm
X Large	WRNTBEXLL	WRNTBEXLR	WRNTBLXLL	WRNTBLXLR	20cm
2X Large	WRNTBE2XLL	WRNTBE2XLR	WRNTBL2XLL	WRNTBL2XLR	21cm

Wrist Braces



KSA ELASTIC WRIST BRACE WITH THUMB

Designed for the control and support of the wrist and thumb. The thumb stay abducts the thumb, supporting the MP joint and restricting movement of the CMC joint.

Size	Left	Right	Wrist
X Small	WRWTXSL	WRWTXSR	14cm
Small	WRWTSL	WRWTSR	16cm
Medium	WRWTML	WRWTMR	18cm
Large	WRWTLL	WRWTLR	19cm
X Large	WRWTXLL	WRWTXLR	20cm



KSA GRADUATED ELASTIC WRIST BRACE

Made from knitted elastic to give a most comfortable feel when worn by shaping to the contours of the wrist and hand. This is our most competitively priced wrist brace. The aluminium stay is covered with polyester webbing. Hook and loop fastening.

Size	Left	Right	Wrist
X Small	WRGRAXSL	WRGRAXSR	14cm
Small	WRGRASL	WRGRASR	16cm
Medium	WRGRAML	WRGRAMR	18cm
Large	WRGRALL	WRGRALR	19cm
X Large	WRGRAXLL	WRGRAXLR	20cm

Note: Also available in sizes 2XL and 3XL



KSA NEOPRENE THUMB SPICA

This neoprene thumb spica is designed to restrict and immobilise movement at the CMC joint and the MCP joint whilst allowing movement at the IP joint and leaving the palm free.

Available in blue and beige.

Size	Blue left	Blue right	Beige left	Beige right	Wrist measurement
Extra small	WRSPIBLXSL	WRSPIBLXSR	WRSPIBEXSL	WRSPIBEXSR	14cm
Small	WRSPIBLSL	WRSPIBLSR	WRSPIBESL	WRSPIBESR	16cm
Medium	WRSPIBML	WRSPIBLMR	WRSPIBEML	WRSPIBEMR	18cm
Large	WRSPIBLLL	WRSPIBLLR	WRSPIBELL	WRSPIBELR	19cm
X Large	WRSPIBLXLL	WRSPIBLXLR	WRSPIBEXLL	WRSPIBEXLR	20cm

Fracture Braces



OVER THE SHOULDER HUMERUS FRACTURE BRACE

For treatment of selected humeral fractures this lightweight fracture brace aims to reduce distal migration and allows full range of movement. Please specify left or right.

Code	Size	Circ	Med Length	Lat Length
RA510552	Med	11"-14.5" (28-37cm)	7" (17.5cm)	12" (30.5cm)
RA510553	Large	12"-16.5" (31-42cm)	7.25" (18.5cm)	12.5" (32cm)
RA510554	XLarge	12"-17.5" (31-42cm)	8" (20.5cm)	13.5" (34.5cm)



MIAMI PREFABRICATED HUMERUS FRACTURE BRACE

Controls humeral segments for ideal treatment of shaft fractures. Still allows the full range of movement in the shoulder.

Code	Size	Circ	Med Length	Lat Length
RA510545SL	Small Left	8-11"(20.5-28cm)	5.25"(13.5cm)	6.25"(16cm)
RA510545SR	Small Right	8-11"(20.5-28cm)	5.25"(13.5cm)	6.25"(16cm)
RA510546ML	Med Left	10-13"(25.5-33cm)	5.5"(14cm)	6.75"(17cm)
RA510546MR	Med Right	10-13"(25.5-33cm)	5.5"(14cm)	6.75"(17cm)
RA510547LL	Lge Left	12-15"(30.5-5.38cm)	6"(15cm)	7.25"(18.5cm)
RA510547LR	Lge Right	12-15"(30.5-5.38cm)	6"(15cm)	7.25"(18.5cm)
RA510548XLL	XLge Left	14-17"(35.5-5.43cm)	6.25"(16cm)	7.75"(19.5cm)
RA510548XLR	XLge Right	14-17"(35.5-5.43cm)	6.25"(16cm)	7.75"(19.5cm)



MIAMI PREFABRICATE TIBIAL FRACTURE BRACE

Soft tissue compression controls tibial segments. Allows full range of movement through the knee and ankle preventing the stiffness that can be associated with immobilization. Measure circumference at largest part of calf.

Code	Description	Size
RA60601SL	Small Left	11.5"-14"(29-35.5cm)
RA60601SR	Small Right	11.5"-14"(29-35.5cm)
RA60602ML	Medium Left	12.5"-15.5"(31.5-39.5cm)
RA60602MR	Medium Right	12.5"-15.5"(31.5-39.5cm)
RA60603LL	Large Left	13"-16.5"(33-42cm)
RA60603LR	Large Right	13"-16.5"(33-42cm)



MIAMI TIBIAL FRACTURE PTB BRACE

Miami Tibial Fracture PTB Brace is a Sarmiento design. Permits joint mobility and early return to normal activities. Optimal control and molding through soft-tissue compression. Comfortable 1/4" polyethelene foam lines anterior shell. Brace suspension and ankle stability controlled with single axis footplate. Water-resistant. Foot component available (sold separately).

		Mid-calf	Patella to Ankle
RA60594L	Small Left	11"-13" (28-33cm)	15" (35.5cm)
RA60594R	Small Right	11"-13" (28-33cm)	15" (35.5cm)
RA60595L	Medium Left	14"-16" (35.5-40.5cm)	16"-17" (40.5-43cm)
RA60595R	Medium Right	14"-16" (35.5-40.5cm)	16"-17" (40.5-43cm)
RA60597L	Large Left	16"-18" (40.5-45.5cm)	17"-18" (43.5-45.5cm)
RA60597R	Large Right	16"-18" (40.5-45.5cm)	17"-18" (43.5-45.5cm)
RA60596L	Left, Med, Long	14"-16" (35.5-40.5cm)	18"-19" (45.5-48cm)
RA60596R	Right, Med, Long	14"-16" (35.5-40.5cm)	18"-19" (45.5-48cm)
RA60598L	Left, Large, Long	16"-18" (40.5-45.5cm)	19"-21" (48-53cm)
RA60598R	Right, Large, Long	16"-18" (40.5-45.5cm)	19"-21" (48-53cm)
RA60604L	Foot Component, Small Left - Womens 4-7, Mens 4-6		
RA60604R	Foot Component, Small Right - Womens 4-7, Mens 4-6		
RA60605L	Foot Component, Medium Left - Womens 7+, Mens 6-9		
RA60605R	Foot Component, Medium Right - Womens 7+, Mens 6-9		
RA60606L	Foot Component, Large Left - Mens 10-12		

Sample Request Form

Name:Department:.....

Hospital:.....

Address:.....

.....

.....

Product Code	Description
1. <input type="text"/>	<input type="text"/>
2. <input type="text"/>	<input type="text"/>
3. <input type="text"/>	<input type="text"/>
4. <input type="text"/>	<input type="text"/>
5. <input type="text"/>	<input type="text"/>
6. <input type="text"/>	<input type="text"/>
7. <input type="text"/>	<input type="text"/>
8. <input type="text"/>	<input type="text"/>
9. <input type="text"/>	<input type="text"/>
10. <input type="text"/>	<input type="text"/>

notes

notes

notes



Kettering Surgical Appliances Limited:
73 Overstone Road
Northampton NN1 3JW.

Telephone: 01604 622886
E-mail: paul@ketteringsurgical.co.uk
Company Registration No. 01019458
Facsimile: 01604 629689
Web: www.ketteringsurgical.co.uk
VAT Registration No. 121 5091 15

KSA